


5 Return

To ship, place all labeled tubes into the foam sleeve. Then place the foam sleeve into the biohazard bag and seal it.


Insert the completed test requisition form and/or the sample info card into the second pocket in the biohazard bag.

Put the sealed biohazard bag into the GeneDx box and close it.

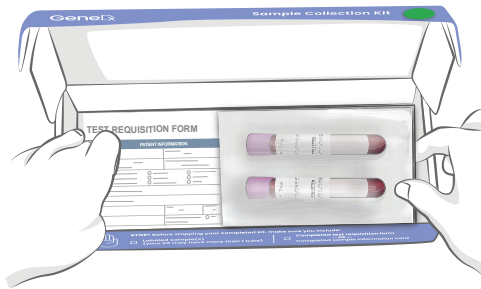

Place the box inside the FedEx Clinical Pak, remove the adhesive strip, and seal the FedEx Pak. Attach the FedEx label to the outside of the Pak.

Call FedEx at (800) 463-3339 to schedule a pickup, or drop off at a local FedEx location. You can find the nearest FedEx location to you by calling FedEx or visiting [FedEx.com/locate](https://www.fedex.com/locate).

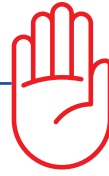

IMPORTANT

Before sealing and shipping the completed kit, make sure to include:



Samples labeled with name, date of birth, and date of collection (2 tubes/individual)

Unlabeled samples will not be processed.



Completed test requisition form

—OR—

Completed sample info card

GeneDx

Blood Sample Collection

Please read these full instructions before collecting the sample.

Take everything out of the GeneDx box and put on a table/counter in front of you.

1 Collect

Perform venipuncture using enclosed tubes. Fill each tube with 2-4 ml of whole blood.



GeneDx


207 Perry Parkway, Gaithersburg, MD 20877
T 1 (888) 729-1206 or 1 (301) 519-2100
F 1 (201) 421-2010

support@genedx.com • www.genedx.com

©2022 GeneDx, LLC. All Rights Reserved



2 Mix

Gently invert each tube 4-5 times to mix the additive with blood.

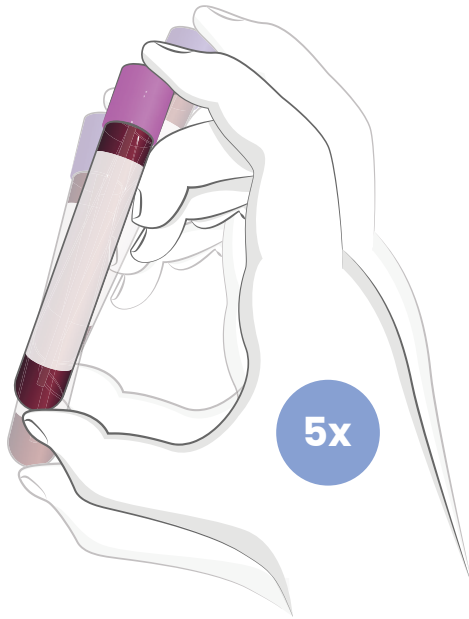

3 Label

Fill out the provided labels with the name and date of birth of the person whose sample is contained in the tubes, as well as the date each sample was collected. Place the labels onto each of the tubes.

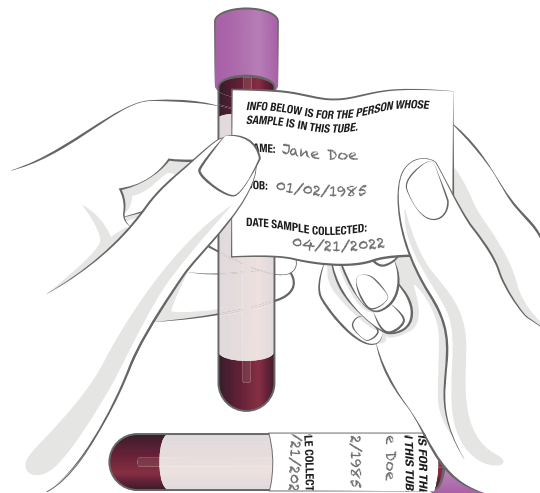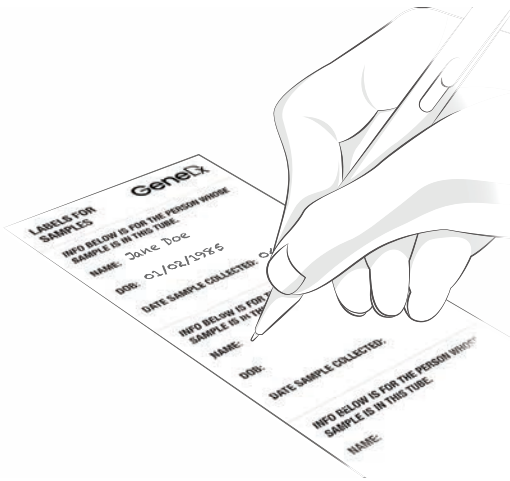


Reminder: Samples submitted without labels will not be processed.

4 Sample Information

Every sample **MUST** be accompanied by a completed test requisition form (preferred) or a completed sample info card.



Test requisition form (preferred)

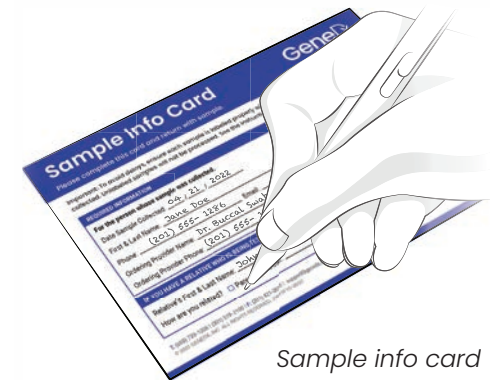

Sample info card

Reminder: Results for samples submitted without this form or card will not be processed.